



RSV Pressure Transmitter



The Anderson RSV Pharmaceutical Series Pressure Transmitter offers not only a new low profile design, but may be ordered with built-in LCD indication. The unit meets all pharmaceutical requirements for performance, finish and cleanability. It is designed and manufactured to withstand the harsh process and sterilization conditions encountered in the pharmaceutical and biotechnology industries.

Its one-piece stainless steel package design incorporates a transducer and electronic circuitry to convert pressure and/or vacuum to a proportional 4-20 mA signal. The wide variety of fittings and ranges provide flexibility in specifying the best transmitter for any application. The unit is designed to operate at high process temperatures and withstand CIP/

SIP conditions.

All units are supplied factory calibrated to standard or custom ranges. For field maintenance, non-interactive zero, and span adjustments, as well as field accessible test points are provided. The optional LCD cap mounted display can be scaled to match actual process units, 0-100% of full scale, or 4-20 mA. This useful feature provides indication directly at the process, in addition to the standard 4-20 mA output. Permanently marked tagging is provided at no extra charge. In addition the RSV complies with Class 1, Div. 1, Groups A, B, C and D for intrinsically safe apparatus. Units are supplied with certificates of calibration.

SPECIFICATIONS

Excitation:

9-32 VDC (Absolute), 24 VDC Nominal regulated or unregulated
14-37 VDC (Absolute) with display

Output:

4-20 mA DC

Loop Resistance:

0-750 ohms at 24 VDC
0-1050 ohms at 30 VDC

Indication:

Optional, 3 1/2 digit, .5" high LCD, cap mounted

Accuracy:

± 0.5% of span, full scale, for standard ranges except those below
0-50 psig (±1.0% full scale. All vacuum/pressure ranges and PSIA ranges are ±1.0%)

Over-Range Rating:

2 times base range

Stability:

± 0.3% of calibrated range / 6 months

Zero and Span Adjustments:

± 10% of range

Repeatability:

± 0.3%

Hysteresis:

± 0.1%

Linearity:

± 0.1%

Process Temperature Range:

20° to 275°F (Up to 325°F if horizontal mount)

Process Temperature Effect:

± 0.1 psig/5.5°C (10°F)

Ambient Temp. Operating Range:

-17.8° to 48.9°C (0° to 120°F)

Storage Temperature:

-40° to 65°C (-40° to 149°F)

Mounting:

Direct connection

Housing Material:

304 Stainless Steel

Housing Ratings:

NEMA 4X, IP66

Wetted Parts:

316L Stainless Steel standard; Hasteloy "C" diaphragm optional

Surface Finish (wetable parts):

R_a max = 8 microinches (.2 microns) electropolished

Recommended Cable:

18-24 AWG, foil shielded, and PVC coated.
(3/16 - 1/4 OD insulation)

Wiring Connection:

Screw Terminal; Accessible via removable screw cap conduit housing

Approvals:

UL Approved; Complies with Class 1, Div. 1, Grps. A-D for intrinsically safe apparatus

Standard:


Designed and manufactured to sound engineering practices in accordance with Article 3.3 of the PED 97/23/EC

Response Time:

200 µSec

* Specifications are subject to change without notice.

Complete Product Ordering Matrix

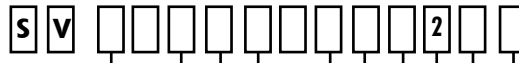

	RANGE		CABLE LENGTH
025	30/0	(M)	00 No Cable
028	30"Hg/0/15 psig	(C)	05 25'
029	30"Hg/0/30 psig	(C)	10 50'
031	30"Hg/0/60 psig	(C)	15 75'
032	30"Hg/0/100psig	(C)	20 100'
066	0/30	(G, A)	21 125'
068	0/50	(G, A, B)	22 150'
069	0/60	(G)	23 175'
071	0/100	(G, A)	24 200'
073	0/150	(G, A)	A0 Quick Disconnect Receptacle(QDR)
074	0/160	(G, A)	A1 QDR & Field Wireable Connector(FWC)
077	0/300	(G, A)	A2 QDR & FWC w/25 ft Cable
081	0/500	(G)	A3 QDR & FWC w/50 ft Cable
084	0/1000	(G)	A4 QDR & FWC w/75 ft Cable
251	-1/0/1 BAR	(B)	A5 QDR & FWC w/100 ft Cable
286	-1/0/2.5 BAR	(B)	A6 QDR & FWC w/125 ft Cable
217	-1/0/3 BAR	(B)	A7 QDR & FWC w/150 ft Cable
056	-1/0/4 BAR	(B)	A8 QDR & FWC w/175 ft Cable
304	-1/0/7 BAR	(B)	A9 QDR & FWC w/200 ft Cable
235	0-3 BAR	(B)	
192	0-4 BAR	(B)	
061	0-10 BAR	(B)	
065	0-20 BAR	(B)	
	ENGINEERING UNITS		
G	PSIG (Gauge)		
M	Hg (inches Mercury - Vacuum)		
A	PSIA (Absolute)		
B	BAR (Gauge)		
C	Vacuum/Pressure Compound (Hg & PSIG)		
T	BAR (Absolute)		
	FITTING		
002	3/4" Tri-Clamp®		
004	1 1/2" Tri-Clamp®		
005	2" Tri-Clamp®		
010	1 1/2" APC "K"		
011	2" APC "K"		
016	1 1/2" Cherry "I" (male)		
017	2" Cherry "I" (male)		
027	1 1/2" G & H "H" line (male)		
028	2" G & H "H" Line (male)		
123	Anderson "CPM" Flushmount		
	DISPLAY VALUE		
0	None		
G	PSIG		
A	PSIA		
B	BAR		
R	mA (milliamphere, for vac/press.)		
Q	% (percent, for vac/press.)		
	DIAPHRAGM MATERIAL		
1	316L Stainless	NOT AVAILABLE WITH FITTINGS "002" or "161"	
2	Hastelloy "C"		
	DIAPHRAGM FINISH		
2	Electropolished (R _a 8)		

SPECIFY TAGGING INFORMATION TO BE PERMANENTLY MARKED, IF APPLICABLE



Indianapolis • Chicago • San Juan
www.hollandapt.com
800-800-8464



ANDERSON INSTRUMENT CO., INC.
156 Auriesville Road
Fultonville, NY 12072 USA
1-800-833-0081 • FAX 516-922-8997
ANDERSON INSTRUMENT CO. LP
400 Britannia Road East, Unit 1
Mississauga, Ontario L4Z1X9 Canada
905-568-1440 • FAX 905-568-1652
NEGELE MESSTECHNIK GmbH
(A Division of Anderson)
Raiffeisenweg 7
D-87743 Egg a. d. Günz Germany
+49 (0) 8333/9204-0 • FAX +49 (0) 8333/9204-49
LATIN AMERICA AND ASIA
Maida International
516-676-3079 • FAX 516-676-3199
www.andinst.com

MEETS OR EXCEEDS CGMP'S OF THE PHARMACEUTICAL INDUSTRY

Form AIC 3562 © 2/04
Revised: December 2008
Supersedes: July 2008